LIMITED WARRANTY

To the consumer, Guardian Technologies LLC warrants this product to be free of defects in materials or workmanship commencing upon the date of the original purchase. Save your original sales receipt to validate start of warranty period. Warranty is not valid without receipt.

If this product should become defective within the warranty period, we will repair or replace any defective parts free of charge. All warranty repairs must be completed by Guardian Technologies LLC. This warranty **does not cover unauthorized repairs**. The warranty does not include unusual wear, damage resulting from accident, or unreasonable use of the product. This warranty only covers the product when used with **genuine Guardian Technologies accessories**. This warranty covers product that was purchased from authorized distributors. This warranty gives you specific legal rights and you may also have other rights (other rights may vary from state to state in the U.S.A.)

The product warranty registration can be completed online at www.guardiantechnologies.com. We consider the registration process important to ensuring superior service to our customers, however submitting the warranty registration is optional and does not affect your rights to utilize this warranty according to the conditions stated above. To submit product under warranty the complete machine must be delivered pre-paid to Guardian Technologies LLC. Please include complete information including: the problem, the model number of the product, the day of purchase, and a copy of the original sales receipt along with your name, address, and telephone (email optional). Address returns to the attention of: Customer Service, at the address below. Additional questions or comments can be made toll free to the number listed below.

**Guardian Technologies LLC
26251 Bluestone Blvd.
Euclid, Ohio 44132
1.866.603.5900
www.guardiantechnologies.com**

Additional copies of the instruction manual can be obtained by accessing and printing directly from our website or calling customer service.


Guardian Technologies LLC
26251 Bluestone Blvd. • Euclid, OH 44132
1.866.603.5900 • www.guardiantechnologies.com

©2016 Guardian Technologies LLC

USE & CARE INSTRUCTIONS

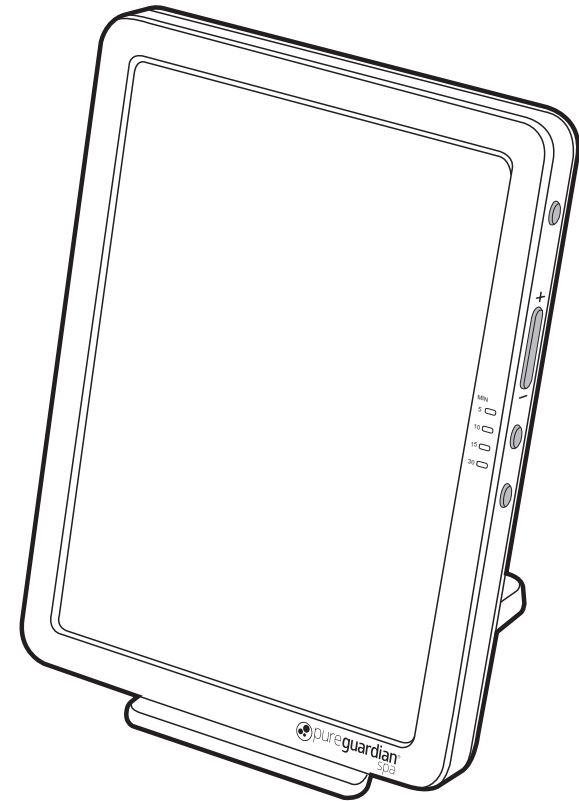Energy Light

Model No. SPA50
1-Year Limited Warranty
www.guardiantechnologies.com
1.866.603.5900
REV1116

English E - 1
French F - 1
Spanish S - 1

Date Purchased
month _____
year _____

READ AND SAVE THESE INSTRUCTIONS

WARNINGS

Failure to comply with the warnings listed below may result in electric shock or serious injury.

This product should be used only in accordance with the specifications outlined in this manual. Usage other than what has been specified here may result in serious injury.

Warnings (a warning indicates the possibility of injury to the user or operator)

- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings (e.g. in the bathroom or near a shower or swimming pool). Do not immerse appliance in water or any other liquid.
- Clean only with a soft, dry cloth – Do not use solvents or cleaners containing abrasives or ammonia.
- Check if the voltage indicated on the appliance corresponds to the local voltage before you connect the appliance.
- The adapter contains a transformer. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- Never use the adapter if it is damaged in any way. Always have the adapter replaced with one of the original type in order to avoid a hazard.
- Never open the appliance to prevent electric shock, finger burns or temporary black spots in the eye.
- If the power cord is damaged, do not use the appliance.
- Close supervision is necessary when this appliance is used by, or near children or invalids.
- Children should be supervised to ensure that they do not play with the appliance.
- There are no user serviceable parts inside the appliance. No modification of this appliance is allowed.
- Never use the appliance if the lamp housing or cover is damaged, broken or missing.
- Do not use the appliance where the light may compromise you or other people's ability to perform essential tasks, such as driving or operating heavy equipment.
- Do not use this appliance in the presence of a flammable anesthetic mixture with air, oxygen or nitrous oxide.
- Consult your doctor before you starting using if:
 - You suffer or have suffered from severe depression;
 - You have an eye condition for which our physician advised you to avoid bright light;
 - You are on specific medication (i.e. certain antidepressants, psychotropic drugs or malaria tablets)

WARNING: This product contains chemicals known to the State of California to cause cancer and birth defects or other reproductive harm. Wash hands after using. (California Code of Regulations Prop 65).

READ AND SAVE THESE INSTRUCTIONS

CAUTION

Caution (a caution indicates the possibility of damage to the appliance)

- Always unplug the appliance after use.
- Do not use the appliance outdoors or near heated surfaces or heat sources.
- Place the appliance on a stable, level and non-slippery surface with the stand in place.
- Do not place any object on top of appliance, adapter or cord.
- When unplugging the appliance, grasp the plug directly to avoid damaging the cord. Never pull on the power cord to remove the plug from the outlet.
- Use this appliance for its intended household use as described in this manual.
- If this appliance is damaged in any way or it is not working it should not be operated.

Do not use this device if:

- You have bipolar disorder
- You are hypersensitive to light and/or experience eye discomfort
- You have had recent eye surgery or a diagnosed eye condition for which your physician advised you to avoid bright light. Or if you have a history of eye disease including, but not limited to, cataracts or macular degeneration.
- You are using any medication that might increase your sensitivity to light (e.g. certain antidepressants, psychotropic drugs or malaria tablets)

If any of the above conditions apply to you, always consult your physician before you start using the device.

Possible Reactions

Possible adverse reactions to light from this device are relatively mild and of a transient nature.

They may include: • Headache • Eyestrain • Nausea

Measures to avoid, minimize or alleviate these adverse reactions:

- Always use the device in a well-lit room.
- Reduce the light intensity if these reactions occur.
- Stop using the device for a few days to make the reactions disappear and then try again.

Intended Use

The Guardian Technologies Energy Light is intended to make people feel more energetic, adjust circadian rhythm, and to provide mood relief from the winter blues.

How Light Influences the Body

Many people suffer from diminished amounts of energy during the winter months, this may be due to the lack of sunlight that you experience. These same side effects can be experienced in normal business environments that have low lighting. Light and its intensity stimulate the receptors in our eyes which can affect how we feel. Bright lights have been used successfully in providing relief for seasonal mood swings.

The Energy Light is an effective alternative to make up for the lack of light that you might experience in day to day life. The device emits bright light to simulate natural sunlight without all of the harmful UV rays.

Duration and Intensity

The National Institute of Health have demonstrated that exposure to white light in the range of 2,500 to 10,000 LUX have helped to reverse the symptoms of "Winter Blues". Exposure time will vary depending on the person. Typical exposure is from 15 minutes to 30 minutes, once or twice a day, depending on individual needs and the intensity and color light selected on the therapy light. The average length of a session of 10,000 lux illumination (white light/high intensity) can be, for example, a much shorter time than for 2,500 lux (blue light/medium intensity). It is recommended to use the device on a regular bases at a consistent time and distance. The user can select between the higher intensity white or lower intensity blue options depending on preference.

When to Use

Depending on when you need the energy boost it might be beneficial to use it at different times. Many people respond best first thing upon awakening. Some, however, do better with evening light. It is necessary to determine the optimum time of day for each individual. It is recommended not to use the light within an hour of going to bed.

PRODUCT SPECIFICATIONS AND PARTS

SPECIFICATIONS

MODEL NUMBER: SPA50



DIMENSIONS: 5.125" L x 6.75" W x 9.5" H (working position)

WEIGHT: Unit - 1.28 lbs. Adaptor - 0.35 lbs

INPUT VOLTAGE: 120V 60Hz 0.2A

USB OUTPUT: DC 5V MAX. 1.5A


LIGHT OUTPUT: White Light 100% Intensity: 10,000 LUX at 7" distance
Blue Light 100% Intensity: 4,500 LUX at 7" distance

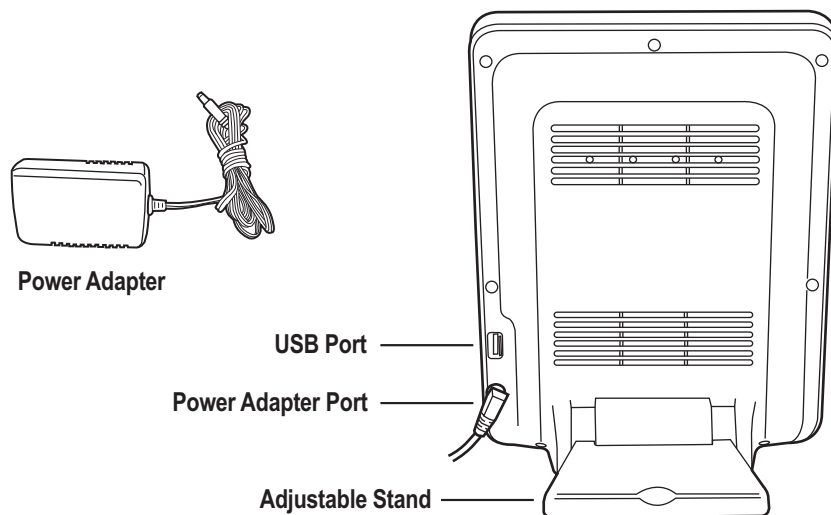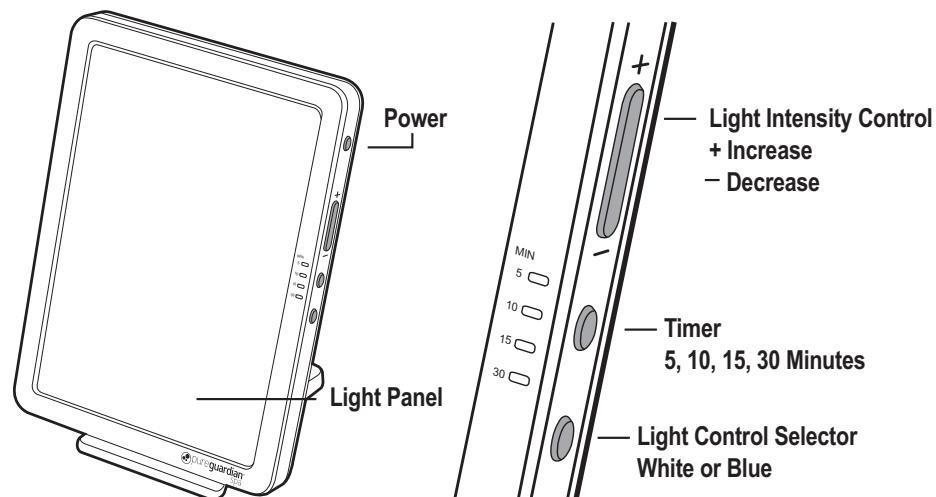INSTRUCTIONS

SET-UP

1. Place the device on a solid surface.
2. Insert the power adapter into the port of the device.
3. Insert the power adapter into an electrical outlet.
4. Placement of the device should be 12-16" away. Adjust the stand at an angle between approximately 90 and 70 degrees.
5. Orient the head and body toward the light. Concentrate on the activities illuminated by the lights (e.g. reading, eating, and computer work) and not on the light itself.
DO NOT STARE DIRECTLY INTO THE LIGHT OF THE DEVICE.
 To minimize eyestrain make sure the room is well-lit.

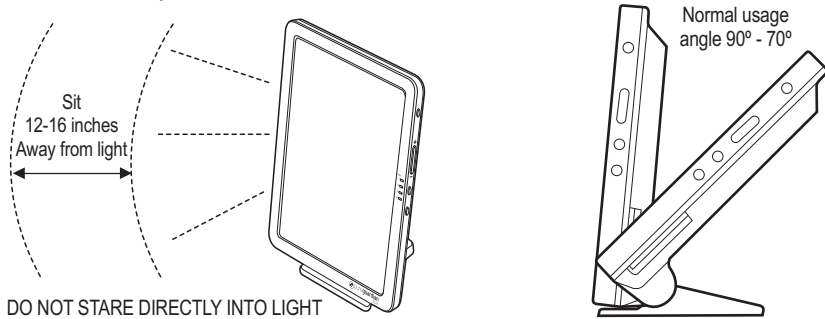USAGE

1. Turn on the device by pressing the power control.
2. Select brightness intensity. There are 3 intensity levels depending on preference.
 - a. Press "+" to increase intensity.
 - b. Press "-" to decrease intensity.
3. Select length of time for session. Push "Timer" control button to toggle between 5, 10, 15 or 30 minutes. Device automatically turns off after selected time. Default time is 5 minutes if no time is selected.
4. Select color. Push "White/Blue" control button and select based on preference. White 2,550-10,000 lux. Blue tint provides range of 1,400-3,750 lux.
5. A USB port is built-in for your convenience to charge other devices while using the energy light.

Cleaning

Clean only with a soft dry-cloth. Never use scouring pads, aggressive liquids, or cleaning agents that contain abrasives or ammonia.

Never immerse the device in water or any other liquid, nor rinse it under the tap.

Storage

If the appliance will not be used for an extended period of time, remove the cord, fold in the base and store in a safe, dry location where it will not be subject to damage.

Disposal

- Dispose of the device in accordance with local regulations.
- Your local or national recycling organizations may also have disposal information.

FAQ's

Question	Answer
Can I get sunburn from using the device?	No, this device does not emit any UV light so it cannot cause sunburn.
Do I have to remove my glasses when I use the device?	No, glasses that are not tinted will not impact the effectiveness of the device.
Can the light harm children or pets?	No, the light is safe to be around children and pets but do not let children play with the device.
When will I start to notice improvement from using the device?	You might experience the energizing effects from the first day or 1 to 2 weeks of daily use.
Can the lamp of the light be replaced?	No, the device uses a durable LED lamp that lasts many years and is not replaceable.

More FAQ's can be found at www.guardiantechnologies.com

TROUBLESHOOTING GUIDE

PROBLEM	SOLUTION
Unit does not operate or shuts off unexpectedly.	<p>Check that the power cord is plugged into the unit.</p> <p>Check that the power cord is plugged into a working electrical outlet.</p> <p>Check that the power button is on.</p> <p>Check that the timer has not turned the unit off.</p> <p>Call Customer Service 1.866.603.5900</p>

FC This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This Class B digital apparatus complies with the Canadian ICES-003.

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